

AllerSnap™ Protein Residue Test



AllerSnap™ Protein Residue Test verifies that proteins have been removed from surfaces. Fast results allow for immediate corrective action to be taken when necessary.

Rating: Not Rated Yet

Price

[Ask a question about this product](#)

Description


AllerSnap™ Protein Residue Test is a ready to use, high sensitivity rapid test for the detection of all protein residues on surfaces. As the system can detect proteins at very low levels, it provides a cost effective tool to screen for the presence of allergens. AllerSnap™ is a user friendly system, providing easy to interpret results within 15 - 30 minutes. A green colour means clean, and a purple colour means re-clean. The more protein present, the quicker the colour change to purple and the darker the colour.



Benefits:


- Results in 15-30 minutes
- Easy-to-interpret colour results
- Broad-range, non specific protein detection offers speed over specific allergen test
- Unique liquid-stable chemistry consistent results
- Pre-moistened swab provides reliable collection, recovery and detection to ensure consistent results

AllerSnap Procedure:

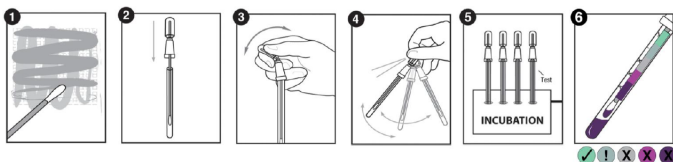

Sensitivity:

- Detects 3µg protein after 30minutes at 37°C or 15 minutes at 55°C.

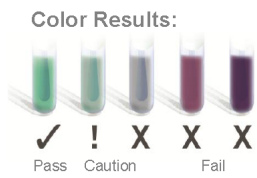

Ordering Details:

- AllerSnap™ Protein Residue Test Product Code: HG-ALS100